

Precision and Accuracy of Axygen® 300 µL Automation Tips for Hamilton® Microlab® STAR™ Line and NIMBUS®



SnAPPShots

A brief technical report from the Corning Applications Group

Srividya Dadi, Audrey Bergeron, and Hannah Gitschier
Corning Incorporated, Life Sciences
Kennebunk, Maine

Introduction

Automated liquid handling and high throughput screening (HTS) are widely used for drug discovery, molecular biology, and genomics. For HTS, reliable sample preparation and delivery methods have become critical to assay performance. Corning's 300 µL Axygen® pipet tips have been specifically designed for applications using the Hamilton Microlab® STAR Line and NIMBUS® liquid handling workstation. The precision and accuracy testing was performed using the Microlab STAR; however, the NIMBUS line of liquid handlers also use the same pipetting system, and, therefore, the tips are compatible with both lines.

The focus of this study was to evaluate the dispensing volume accuracy and precision of the Axygen 300 µL tips on the Hamilton Microlab STAR automation platform as compared to Competitor 300 µL tips. These criteria were measured using the Artel Multichannel Verification System (MVS®), which calculates the volume of dispensed samples with an absorbance-based measurement system. The results demonstrate that Axygen 300 µL tips are comparable to Competitor 300 µL tips using the Hamilton Microlab STAR liquid handling workstation to dispense volumes as low as 30 µL and as high as 300 µL.

Materials/Methods

Materials

Axygen 300 µL tips (Corning Cat. No. HT-300-CBK-HTR) and Competitor 300 µL tips.

Methods

The Hamilton Microlab STAR liquid handling workstation (Hamilton Cat. No. 1532) was used to assess accuracy as percent deviation (% D) and precision as coefficient of variation (% CV), for Axygen 300 µL tips and Competitor 300 µL tips.

To test the ability of each brand of tips to dispense accurately and precisely, a column of 8 tips was arranged so that each tip aspirated from an Axygen low profile reservoir (Corning Cat. No. RES-SW96-LP) and dispensed into 1 column of a Corning 96-well black clear-bottom microplate (Corning Cat. No. 3631).

For the 30 µL test volume, each tip aspirated 30 µL of Range B solution (Artel Cat. No. MVS-204) and dispensed 30 µL into 170 µL of diluent solution (Artel Cat. No. MVS-202) in each well. For the 300 µL test volume, each tip aspirated 300 µL of Range HV solu-

tion (Artel Cat. No. MVS-214) and dispensed 300 µL into each well. To determine the volume of liquid dispensed into each well, absorbance readings for the solutions (diluted Range B solution for 30 µL dispense and Range HV solution for 300 µL dispense) were measured using an Artel ELx800NB® Plate Reader (Artel Cat. No. 1311197). Each study was performed 6 independent times for each brand of tips for a total of 48 tip dispenses. Evaluation criteria include % D from the set dispense volume and the % CV for the 48 tip dispenses.

Results/Discussion

The evaluation criteria for comparing Axygen 300 µL tips with Competitor 300 µL tips are listed in Tables 1 and 2. The ability of the pipet tips to dispense 30 µL and 300 µL volumes accurately and precisely was determined through the analysis of the mean volume dispensed across 48 wells. The precision of each brand of tip is represented by the % CV of the replicates. Similarly, the accuracy is represented by the % D from the target volume of the replicates. It is important to note that the accuracy of liquid dispense may vary depending on the method and liquid class selection chosen when using the liquid handling platform. However, for these studies the method and liquid used for testing was identical for Axygen 300 µL tips and Competitor 300 µL tips.

Table 1. Evaluation Criteria for 30 µL Dispense Volume

30 µL	Axygen	Competitor
n	48	48
Target Volume (µL)	30.00	30.00
% CV	0.39% ± 0.08%	0.49% ± 0.15%
% D	3.43% ± 0.08%	3.25% ± 0.21%
Total No. of Outliers	0	0

Table 2. Evaluation Criteria for 300 µL Dispense Volume

300 µL	Axygen	Competitor
n	48	48
Target Volume (µL)	300.00	300.00
% CV	0.53% ± 0.10%	0.52% ± 0.12%
% D	2.04% ± 0.39%	2.14% ± 0.17%
Total No. of Outliers	0	0

Data in tables shows ± standard deviation

As demonstrated in Figure 1, Axygen® 300 µL tips displayed comparable precision to Competitor 300 µL tips using the Hamilton® Microlab® STAR™ automation system. There was no significant difference in the precision of each brand of tips when dispensing 30 µL (Figure 1A) or 300 µL (Figure 1B).

As demonstrated in Figure 2, Axygen 300 µL tips displayed comparable accuracy to Competitor 300 µL tips using the Hamilton Microlab STAR automation system. There was no significant

difference in the accuracy of each brand of tips when dispensing 30 µL (Figure 2A) or 300 µL (Figure 2B).

Conclusions

Axygen 300 µL tips demonstrate precision and accuracy comparable to competitor 300 µL tips using the Hamilton Microlab STAR liquid handling workstation to dispense volumes as low as 30 µL and as high as 300 µL.

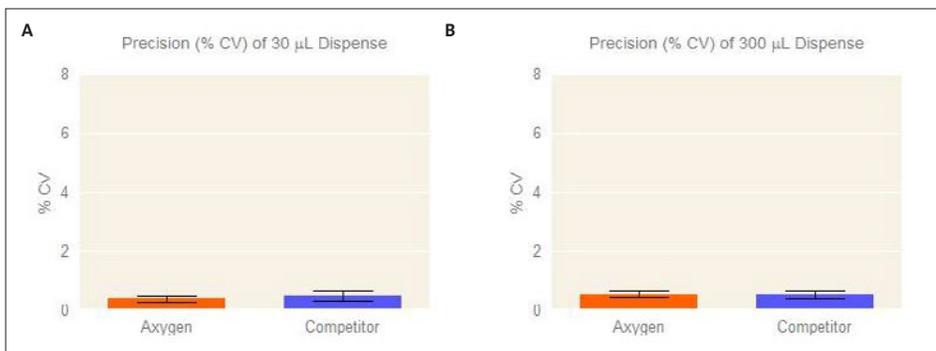


Figure 1. Precision (% CV) analysis of 300 µL tips. The % CV of Axygen and Competitor 300 µL tips dispensing (A) 30 µL and (B) 300 µL volumes using the Hamilton Microlab STAR liquid handling workstation was determined using the Artel MVS® System. There was no significant difference in the % CV between each brand. Data shown with standard deviation (SD). n=48.

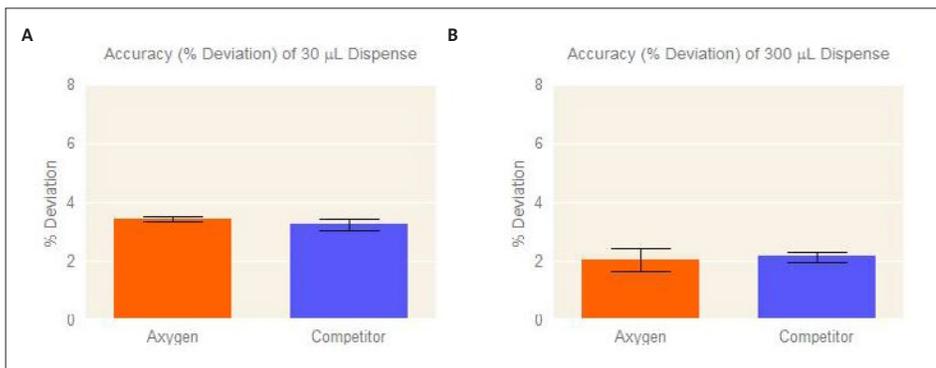


Figure 2. Accuracy (% D) analysis of 300 µL tips. The % D of Axygen and Competitor 300 µL tips dispensing (A) 30 µL and (B) 300 µL volumes using the Hamilton Microlab STAR liquid handling workstation was determined using the Artel MVS System. There was no significant difference in the % D between each brand. Data shown with SD. n=48.

Warranty/Disclaimer: Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications.

For additional product or technical information, please call 800.492.1110 or visit www.corning.com/lifesciences. Customers outside the United States, call +1.978.442.2200 or contact your local Corning sales office listed below.

Corning Incorporated Life Sciences

836 North St.
Building 300, Suite 3401
Tewksbury, MA 01876
t 800.492.1110
t 978.442.2200
f 978.442.2476

www.corning.com/lifesciences

Worldwide Support Offices

ASIA/PACIFIC
Australia/New Zealand
t 0402-794-347

China
t 86 21 2215 2888
f 86 21 6215 2988

India
t 91 124 4604000
f 91 124 4604099

Japan
t 81 3-3586 1996
f 81 3-3586 1291

Korea
t 82 2-796-9500
f 82 2-796-9300

Singapore
t 65 6733-6511
f 65 6861-2913

Taiwan
t 886 2-2716-0338
f 886 2-2516-7500

EUROPE

France
t 0800 916 882
f 0800 918 636

Germany
t 0800 101 1153
f 0800 101 2427

The Netherlands
t 31 20 655 79 28
f 31 20 659 76 73

United Kingdom
t 0800 376 8660
f 0800 279 1117

All Other European Countries

t 31 (0) 20 659 60 51
f 31 (0) 20 659 76 73

LATIN AMERICA

Brasil
t (55-11) 3089-7419
f (55-11) 3167-0700

Mexico
t (52-81) 8158-8400
f (52-81) 8313-8589

CORNING | **FALCON** | **AXYGEN** | **GOSSSELIN** | **PYREX**

For a listing of trademarks, visit us at www.corning.com/lifesciences/trademarks.

Hamilton, MICROLAB, and NIMBUS are registered trademarks, and CO-RE, MICROLAB STAR, STARlet, and STARplus are trademarks of Hamilton Company. The compatibility of these tip products has been independently determined by Corning and there should be no implication of sponsorship or endorsement by Hamilton Company. All other trademarks in this document are the property of their respective owners.