

# Performance Testing for Axygen® Automation Tip (TTF-1000-CBK-HTR-S)

## Application Note



### Method


The liquid handling arm (LiHa) on the Tecan® Freedom EVO® liquid handling workstation was used to assess precision as coefficient of variation (% CV), and accuracy as percent deviation (% D), for Axygen 1,000 µL tips.

To test the ability of the tip to dispense accurately and precisely at two dispense volumes, 100 µL and 1,000 µL, a column of 8 tips aspirated from an Axygen low profile reservoir (Corning Cat. No. RES-SW96-LP) and dispensed into a Corning® 96-well, black, clear bottom microplate (Corning Cat. No. 3631).

For the 100 µL test volume on the LiHa, a column of 8 tips were arranged so that each tip aspirated 100 µL of Range A solution (Artel Cat. No. MVS-203) and dispensed 100 µL into 100 µL of diluent solution (Artel Cat. No. MVS-202) in each well.

For the 1,000 µL test volume on the LiHa, a column of 7 tips were arranged so that each tip aspirated 1,000 µL of Range HV solution (Artel Cat. No. MVS-214) and dispensed 250 µL into 0 µL of diluent solution in each well. To determine the volume of liquid dispensed in each well, absorbance readings for the solutions (diluted Range A solution for 100 µL dispense and Range HV solution for 1,000 µL dispense) were measured using an Artel ELx800NB® plate reader (Artel Cat. No. 1311197). Each study was performed 6 independent times for a total of 48 tip dispenses for the 100 µL test volume and 3 independent times for a total of 21 tip dispense for the 1,000 µL test volume. Evaluation criteria include % D from the set dispense volume and % CV of the measured dispense volume for the tip dispenses.

### Results

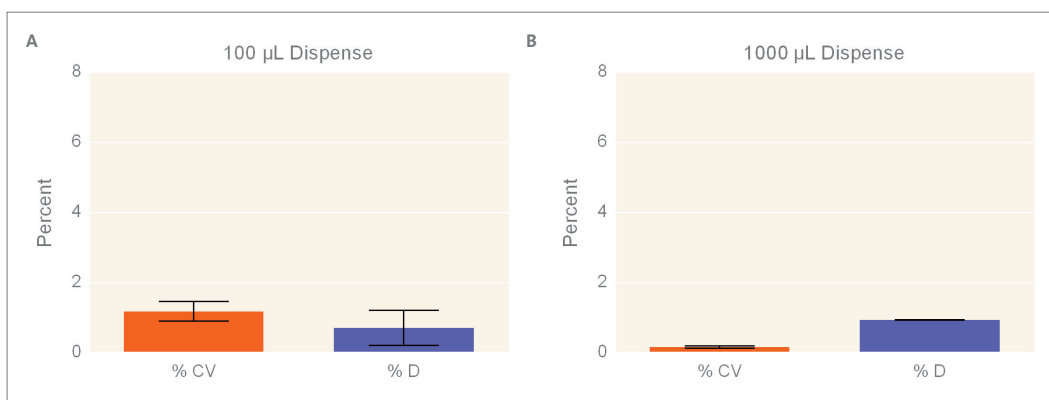


**Figure 1.** Analysis of TTF-1000-CBK-HTR-S tip with aqueous dispense. The precision (assessed by % CV) and accuracy (assessed by % D) of Axygen TTF-1000-CBK-HTR-S tips dispensing (A) 100 µL and (B) 1,000 µL volumes using the LiHa on the Tecan Freedom EVO liquid handling workstation were determined using the Artel MVS® system. The % CV and % D were below 1.5% for both 100 µL and 1,000 µL dispenses.

**Table 1.** Aqueous Dispense Results

Target Volume (µL)	100	1,000
n	48	21
% CV	1.19 ± 0.28	0.17 ± 0.03
% D	0.72 ± 0.50	0.95 ± 0.01
Outliers	0	0

### Conclusion

The % CV and % D for the Axygen automation TTF-1000-CBK-HTR-S tips dispensing 100 µL and 1,000 µL were 5% or below. Therefore, Axygen automation TTF-1000-CBK-HTR-S tips can precisely and accurately dispense volumes as low as 100 µL and as high as 1,000 µL for aqueous solution using the LiHa on the Tecan Freedom EVO liquid handling workstation.

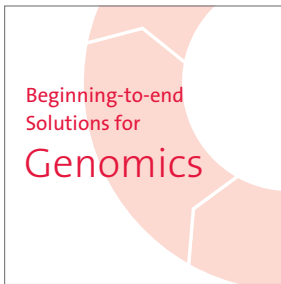

[www.corning.com/lifesciences/solutions](http://www.corning.com/lifesciences/solutions)

In our continuous efforts to improve efficiencies and develop new tools and technologies for life science researchers, we have scientists working in Corning R&D labs doing what you do every day, across the globe. From collection to analysis, our technical experts understand your challenges and your need for simplified efficient, low- to high-throughput genomics processes.

A combination of global manufacturing expertise, extensive use of in-house automation, an unsurpassed commitment to product innovation and a thorough understanding of your processes enables Corning to offer a beginning-to-end portfolio of high-quality, reliable consumables and reagents for genomics applications.

For more specific information on claims, visit the Certificates page at [www.corning.com/lifesciences](http://www.corning.com/lifesciences).  
**Warranty/Disclaimer:** Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures. Not for use in humans. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications.

For additional product or technical information, visit [www.corning.com/lifesciences](http://www.corning.com/lifesciences), or contact our Scientific Support Team at [ScientificSupportEMEA@corning.com](mailto:ScientificSupportEMEA@corning.com).

**Corning Incorporated**  
*Life Sciences Europe*

Corning BV  
Fogostraat 12  
1060 LJ Amsterdam  
The Netherlands  
Phone: +31 (0) 20 659 60 51  
Fax: +31 (0) 20 659 76 73  
[CSEurope@corning.com](mailto:CSEurope@corning.com)  
[www.corning.com/lifesciences](http://www.corning.com/lifesciences)

**Support Offices**

**EUROPE**

**France**  
t 0800 916 882  
f 0800 918 636  
**Germany**  
t 0800 101 1153  
f 0800 101 2427

**The Netherlands**  
t 31 20 655 79 28  
f 31 20 659 76 73  
**United Kingdom**  
t 0800 376 8660  
f 0800 279 1117

**All Other European Countries**  
t 31 (0) 20 659 60 51  
f 31 (0) 20 659 76 73

**CORNING** | **FALCON** | **AXYGEN** | **GOSSELIN**

For a listing of trademarks, visit [www.corning.com/clstrademarks](http://www.corning.com/clstrademarks).  
All other trademarks are the property of their respective owners.