

The Corning® Guide to Surface Selection by Cell Type

The right surface — right from the start

The development and normal functioning of cells depends on interactions with molecules in their microenvironment. To create physiologically relevant *in vitro* models that support normal cell growth and function, the components of the *in vivo* environment must be incorporated. Use of extracellular matrix (ECM) proteins (natural, synthetic, or mimetic) as coatings on the cultureware allows the development of cell type specific model systems which more closely mimic *in vivo* conditions. Use the Corning selection guide to help you choose the best surface for your application and cell type.

The numbers in the following tables refer to the applicable document references, which are listed at the end of this document.

For product information, please refer to the Corning Surfaces Brochure (CLS-C-DL-AC-006) or visit www.corning.com/lifesciences/advancedsurfaces.

CORNING

Primary Cells

Primary Cells	Extracellular Matrices (ECMs) and Biological Coatings													ECM Mimetics and Advanced Surfaces				Enhanced TC-treated Surfaces				
	Cell-Tak™	Collagen I	Collagen IV	Matrigel® Matrix	Fibronectin	Gelatin	Laminin	Osteopontin	Poly-Lysine (PDL, PLL)	PDL/LM and PLO/LM	PuraMatrix®	Vitronectin	PureCoat™ ECM Mimetic Fn	PureCoat ECM Mimetic COL I	Synthemax® Surface	Ultra-Low Attachment	Osteo Assay Surface	laminin-521 (Human)	Primaria™	CellBIND® Surface	PureCoat Amine	PureCoat Carboxyl
Aortic endothelial cells, BAEC		1		2, 3	4		4					5										
Bile duct cells (epithelial)		6		7																		
Bone marrow cells (bone resorption, osteoclast)																8-10						
Brain microvessel (endothelial)		11, 12	12	13	12	14	12					12										
Cardiomyocytes; cardiac (endothelium, progenitor cells)		15		16	17		18		19		20								21			22
Colonocytes (epithelial)			23	24											25							
Dorsal root ganglia				26, 27					28	29												
Embryonic cortical neurons				30							31											
Embryonic sympathetic neurons			32	33			32				32											
Endothelial cells; endothelial colony forming cells			34		34, 35		34						36	36					37			
Erythrocyte culture (parasite development stages [asexual, sexual])	378			379																		
Hepatocytes		38, 386, 387	39, 387	40, 387			387		41		42								43, 44	45		
Hippocampal neurons				46, 47	48		46, 48		46, 49, 50	49, 50	47, 51											
Human periodontium (periodontal ligament)	52																					
Human osteoclast precursors (osteoclast, pit formation)																53, 54						
HUVEC (endothelial)		55, 56		57, 58	56, 59-61	62	56	63			64	56, 58, 61							61			
HVSMC				58			65					58										
Keratinocytes		66, 67		67	67, 68						51	68		315		69						
Mammary epithelial cells; breast cells (luminal, myoepithelial and endothelial)		70, 71, 73		71, 72			74				71					75						
Microvascular, BME (endothelial)		76	77	78	79, 80	76, 81					82	80										
Mouse splenic T-cells	83		84	84																		
Muscle cells, myoblasts, myogenic cells, myotubes				85			86													87		
Neuronal cells (cortical, cerebellar granule, astrocytes, sensory, sympathetic)			88				88, 89		90-98	99											100	
Oligodendrocytes (glial; precursors)				101			102		101, 103			102										
Osteoblasts		104									105	104										
Pancreatic islet, neonatal (3- to 5-day-old) rat islets of langerhans	106			107	106											108						22
Parotid acinar cells	109			110																		
Peripheral blood mononuclear cells		111, 112	113	114	112-114							113				115	116					
Postnatal mouse vestibular ganglion neurons	117																					
Schwann cells (glial)			118	118			118				119											
Sertoli cells (spermogenic)	120			121, 122																		
Skeletal muscle cells (myocytes, myotubes)				123															124	125		
Smooth muscle cells (endothelial, aortic, vascular)	373	126	126	127	126, 128														129			
Urothelial cells		130	130	131	132																	
Valvular interstitial cells					133																	
Zygote and blastocyst development stages	375																					

Cell Lines (transformed or transfected)

Cell Lines	Extracellular Matrices (ECMs) and Biological Coatings											ECM Mimetics and Advanced Surfaces					Enhanced TC-treated Surfaces					
	Cell-Tak™	Collagen I	Collagen IV	Matrigel® Matrix	Fibronectin	Gelatin	Laminin	Osteopontin	Poly-Lysine (PDL, PLL)	PDL/LM and PLO/LM	PuraMatrix®	Vitronectin	PureCoat™ ECM Mimetic Fn	PureCoat ECM Mimetic COL I	Synthemax® Surface	Ultra-Low Attachment	Osteo Assay Surface	rLaminin-521 (Human)	Primaria™	CellBIND® Surface	PureCoat Amine	PureCoat Carboxyl
ARH-77 (lymphoblast)					255																	
BHK-21 (fibroblast)					61	256					61								61		257, 258	
Breast cancer cells (established cell lines)	259, 260			261, 262				259														
C2C12 (myoblast)		263		264							265					266						
Cell immobilization (Gin-1, Nasal epithelial cells, Molt-4 and K562 human leukemia cells, Sf9 Cells)	267-270																					
Chinook Salmon Embryo Cells (CHSE-214)																				271		
CHO, CHO-1, CHO-K1 (epithelial, endothelial, transfected fusion protein)				272, 273				274		51		275						276	277, 278	22		
COS-7 (fibroblast, transfected)		279		280	279			281, 282			279							283				
Dorsal Root Ganglia (transfected)				284							285											
H1299 (transfected- human non-small cell lung carcinoma cell line)				286	287																	
HEK-293 (transfected, epithelial), EcoPack2™-293, HEK-SRATet cells, Living Colors HEK-ZsGreen proteasome sensor (transfected)	288	289, 290		291		291		274, 292		51					293			294	295, 296	258, 297	297	
HeLa										51												22
HepG2 (hepatocyte), Hep3B (hepatoma)		299		300						51	377				301				125	258	258	
HT-1080 (epithelial)		302, 303	304, 305	302											307, 308							258
hFOB 1.19, MG63 (osteoblast cell lines)				309-311	312		313			311	312				314							
Human MOLT-4, drosophila S2 (biomaterial and tissue engineering applications)	374																					
Keratinocytes (human neonatal)		315, 316			316								315									
L929 (fibroblast, transfected)				317			318			319												
LnCAP (prostate cancer cell line)		307		320											307				296		257, 258	
MCF7 (epithelial)		321	322		323						322				324							
MCF-10A (epithelial)		71, 325		71, 325 - 328	329, 330		330, 331		332	71	330				333							
MDA-MB-231		302, 307, 334	322	302, 326, 335-339	322, 334	322	334	334			322				307, 324							
MDA-MB 435		340		338, 339, 341, 342							343											
MM41 (skeletal myoblasts, transfected)		344																				
MRC5																						257
N2AB-1 (neuroblastoma)	385																					
NIH/3T3, 3T3 (fibroblast)				345	346, 347			349		51												
PC-3, PC-12		307, 350		351, 390			352	348, 353, 354	355	47, 51								356	357	22	22	
RTG-2 (rainbow trout gonad cells)				358																271		
RAW 264.7 (macrophage; osteoclast differentiation, pit formation)			359				359									360, 361						
SH-SY5Y		362	362	363			364		364	51	365											
SK-MEL-28			366		366, 367		366				367											
U266 (lymphoblast)					255																	
U937 (monocyte)		368					369				370				371							
Vero cells												275	275									

Stem and Progenitor Cell Expansion

	Cell-Tak™	Collagen I	Collagen IV	Matrigel® Matrix	Fibronectin	Gelatin	Laminin	Osteopontin	Poly-Lysine (PDL, PLL)	PDL/LM and PLO/LM	PuraMatrix®	Vitronectin	PureCoat™ ECM Mimetic Fn	PureCoat ECM Mimetic COL I	Synthemax® Surface	Ultra-Low Attachment	Osteo Assay Surface	rLaminin-521 (Human)	Primaria™	CellBIND® Surface	PureCoat-Amine	PureCoat Carboxyl
Stem and Progenitor Cells	Extracellular Matrices (ECMs) and Biological Coatings												ECM Mimetics and Advanced Surfaces				Enhanced TC-treated Surfaces					
Human embryonic stem cell (hESC)			134	135	134		134					134			135	136		392		134, 137		
Human induced pluripotent stem cell (hiPSC)				138, 139											138			393				
hMSCs (bone marrow derived, adipose derived)					140			140				140	141	140, 142						397		
Human retinal progenitor cells (RPE)					143									143								
rESC; rat endothelial progenitor cells						144						145, 146				144						
Neuronal stem cell (intestinal/enteric)					147		147									147						

In Vitro Differentiation of Pluripotent Stem Cells

Stem Cells	Extracellular Matrices (ECMs) and Biological Coatings												ECM Mimetics and Advanced Surfaces				Enhanced TC-treated Surfaces					
hESC (cerebral organoid model)				148																		
hESC (pancreatic)				149		150																
hESC, hiPSC (cardiomyocytes)				149, 151-153		154									135			394				
hESC, hiPSC, mESC (Germ Cell Layers: ectoderm, mesoderm, endoderm; hematopoietic progenitor; definitive differentiation; cardiomyocytes)		155, 228	156	135, 138, 139, 158, 159	157	160		161, 162				157		135, 138, 159, 163	157			393				
hESC, hiPSC, mESC, miPSC (endothelial)	164	165		164, 166, 167				164														
hESC, hiPSC (intestinal organoids)				168, 169												147						
hESC, hiPSC (neuronal)				149, 170, 171	170		170, 172		170	173	175	170				174, 176		393				
hESC (osteogenic)						177																
hESC, hiPSC (smooth muscle)				170, 178	170		170		170			170										
hESC, mESC (lung epithelial)		179		179, 180		182										181						
hESC, mESC, rESC (hepatocyte, hepatocyte-like)		183		170, 183-186	170	187	170		170			170				183						
Human NPCs (differentiation to neuronal cells)				188			189					188-190						393				
hPSCs, mPSCs (renal progenitor cells, renal tubular cells, endoderm)		155		191, 192												192						
mESC (hematopoietic)	164			164			164															
mESC, Chicken (cardiomyocytes)		165, 193		193	165	154, 194, 195	165															
mESC, rESC, miPSC (neuronal, progenitor)				188	188, 196	197, 198	188, 199		200		188					198						
mPSCs (inner ear sensory epithelia)				201																		
hESC, hiPSC (retinal pigment epithelial)				396											395							

In Vitro Differentiation of Adult Stem Cells

Stem Cells	Extracellular Matrices (ECMs) and Biological Coatings													ECM Mimetics and Advanced Surfaces				Enhanced TC-treated Surfaces				
	Cell-Tak™	Collagen I	Collagen IV	Matrigel® Matrix	Fibronectin	Gelatin	Laminin	Osteopontin	Poly-Lysine (PDL, PLL)	PDL/LM and PLO/LM	PuraMatrix®	Vitronectin	PureCoat™ ECM Mimetic Fh	PureCoat ECM Mimetic COL I	Synthemax® Surface	Ultra-Low Attachment	Osteo Assay Surface	rLaminin-521 (Human)	Primaria™	CellBIND® Surface	PureCoat Amine	PureCoat Carboxyl
hADSCs; adipose (endothelial)				202													203					
Cardiac progenitor cells (cardiomyocyte)		204					205		204			205					206					
Colon (epithelial organoids)		207		169, 208													209					
Hair follicle (melanocytes, neurons, smooth muscle)				210	210																	
Hepatic progenitor cells (hepatic, biliary cells)							211										212					
Intestinal (organoids, crypt-villus)		213		214-216																		
Keratinocytes (epidermal)		217				217																
Lung (sphere)				218													219					
Mammary epithelial cells				220 - 222													222					
MSC (cardiomyocyte, chondrocyte, hematopoietic, hepatocyte, neuron, osteocyte, spheroid)		141, 223-227, 232		223, 225, 229, 230	140, 223-225, 229, 231		223, 224	140				233-237	140, 232				238, 239					
MSC (endothelial progenitors)		240										240					144					
Muscle (skeletal)							241															
Neural progenitor/stem cells (neuron, astrocytes, neuroblast)				243		198	189, 243					244	189, 190, 242, 245-247				198					
Pancreatic (endocrine)			249	248, 250			249															
Prenatal rat cells (neuron, glial cells)							251															
Retinal (retinal neuron)												252										
Salivary gland				253																		
Stomach (gastric units)				254																		

3D Cell Culture Applications

Cells	Extracellular Matrices (ECMs) and Biological Coatings											ECM Mimetics and Advanced Surfaces				Enhanced TC-treated Surfaces							
	Cell-Tak™	Collagen I	Collagen IV	Matrigel® Matrix	Fibronectin	Gelatin	Laminin	Osteopontin	Poly-Lysine (PDL, PLL)	PDL/LM and PLO/LM	PuraMatrix®	Vitronectin	PureCoat™ ECM Mimetic Fn	PureCoat ECM Mimetic COL I	Synthemax® Surface	Ultra-Low Attachment	Osteo Assay Surface	rLaminin-521 (Human)	Primaria™	CellBIND® Surface	PureCoat Amine	PureCoat Carboxyl	
4T1 (mouse breast cancer cell line)				372																			
Cardiac fibroblast		376																					
Hep3B (hepatoma; toxicity/drug screening)		321																					
MCF-7 (epithelial)		321													324								
MCF-10A (epithelial)		71		71, 326						71					333								
MDA-MB-231		302, 307		302, 326											307, 324								
MDA-MB-361				326																			
HeLa				391											298								
HT-1080 (epithelial)		302, 307		302, 303											307								
hESC, Rat (endothelium)		145, 380		146						240					381								
Human melanoma cell lines SBCL2 (RGP), WM-115, (VGP) and 451-LU (MM) and keratinocytes (spheroid model)		382																					
Mouse embryonic pancreatic progenitors (organoid)				383																			
MSCs, Ovarian cancer cells (OCC)				384											384								
Primary rat hepatocytes				387						42													
Rat hepatocyte progenitor cells (spheroid)										42, 388													
SK-MEL-28 cells				306																			
MEFs (stromal fibroblast)				372																			
U266 (lymphoblast)				389																			

The data in this surface selection guide has been derived from published papers accessed through NCBI databases, as well as various web references. This guide will be periodically updated as additional literature becomes available.

Warranty/Disclaimer: Unless otherwise specified, all products are for research use only. Not intended for use in diagnostic or therapeutic procedures. Not for use in humans. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications.

References

1. Koyama S, et al. (2002) *Br J Pharmacol.* 137:901-9.
2. Adamcic U, et al. (2012) *PLoS One.* 7:e53385.
3. Suh SH, et al. (1999) *Pflugers Arch.* 438:612-20.
4. Shimizu H, et al. (1999) *J Biol Chem.* 274:11995-2000.
5. Bernardis MT, and Jiang S. (2008) *J Biomed Mater Res A.* 87:505-14.
6. Auth MK, et al. (2001) *Hepatology.* 33:519-29.
7. Mathis GA, and Siric AE. (1990) *In Vitro Cell Dev Biol.* 26:113-8.
8. DeMambro VE, et al. (2012) *J Bone Miner Res.* 27:390-400.
9. Ha H, et al. (2013) *BMC Complement Altern Med.* 13:112.
10. Zhong Z, et al. (2012) *Proc Natl Acad Sci USA.* 109:E2197-204.
11. Kanda T, et al. (1994) *J. Cell Biol.* 126:235-46.
12. Wang J, and Milner R. (2006) *J Neurochem.* 96:148-59.
13. Doron DA, et al. (1991) *In Vitro Cell Dev Biol.* 27A:689-97.
14. Roberts TK, et al. (2012) *Lab Invest.* 92:1213-33.
15. Bjorkegren J, et al. (2001) *J Biol Chem.* 276:38511-17.
16. Zhou J, et al. (2010) *J Cell Mol Med.* 14:2641-5.
17. Hazeltine LB, et al. (2012) *Int J Cell Biol.* 2012:508294. doi: 10.1155/2012/508294
18. Simpson DG, et al. (1993) *J Cell Biol.* 123:323-36.
19. Sander V, et al. (2013) *Nat Protoc.* 8:800-9.
20. Tokunaga M, et al. (2010) *J Mol Cell Cardiol.* 49:972-83.
21. Deschepper CF, et al. (2001) *Circ Res.* 88:223-8.
22. Corning PureCoat Cultureware. (2009) Corning Brochure. CLS-DL-CC-046 REV1.
23. Desai GN, et al. (1996) *Gastroenterology.* 111:1541-50.
24. Reddy PM, et al. (1996) *Pediatr Res.* 39:287-94.
25. Chen HJ, et al. (2012) *J Clin Invest.* 122:3184-96.
26. Khan Z, et al. (2002) *Neuroscience.* 114:601-9.
27. Liazoghli D, et al. (2012) *ACS Chem Neurosci.* 3:90-95.
28. Moore K, et al. (2006) *Tissue Eng.* 12:267-78.
29. Nakashima K, et al. (1999) *J Neurosci.* 19:5429-34.
30. Kim YG, and Lee YI. (2012) *Exp Neurobiol.* 21:61-7.
31. Liu G, et al. (2004) *Nat Neurosci.* 7:1222-32.
32. Lein PJ, et al. (1991) *J Cell Biol.* 113:417-28.
33. Ghogha A, et al. (2012) *J Vis Exp.* Mar 21, (61): 10.3791/3546 3546
34. Ando J, et al. (1991) *In Vitro Cell Dev Biol.* 27A:525-32.
35. Ingber DE. (1990) *Proc. Natl. Acad. Sci. USA.* 87:3579-83.
36. Partridge J, et al. (2012) Application Note 492. Corning Literature CLS-DL-CC-007.
37. Yue BY, et al. (1989) *Invest Ophthalmol Vis Sci.* 30:248-53.
38. Bajt ML, et al. (2004) *Toxicol Sci.* 80:343-9.
39. Farkas MH, et al. (2003) *J Biol Chem.* 278:9412-7.
40. Harris AJ, et al. (2004) *Mutat Res.* 549:79-99.
41. Sakai Y, et al. (1991) *Biotechnology Techniques.* 6:527-532.
42. Wang S, et al. (2008) *Tissue Eng Part A.* 14:227-36.
43. Barrera A, et al. (2005) *J Virol.* 79:9786-98.
44. Francavilla A, et al. (1986) *Cancer Res.* 46:1318-23.
45. Schacke M, et al. (2009) *J Virol Methods.* 157:25-31.
46. Sun Y, et al. (2012) *Biointerphases.* 7:29.
47. Holmes TC, et al. (2000) *Proc Natl Acad Sci USA.* 97:6728-33.
48. Chen WH, et al. (2013) *PLoS One.* 8:e83394.
49. Maiese K, et al. (1993) *J. Neurosci.* 13:3034-40.
50. Craig AM, et al. (1993) *Neuron.* 10:1055-68.
51. Zhang S, et al. (1995) *Biomaterials.* 16:1385-93.
52. Shimizu Y, et al. (1993) *J Oral Implantol.* 19:14-22.
53. Faruqi AF, and Beltzer J. (2010). Corning Literature CLS-AN-162.
54. Hongwei, R, et al. (2010). Corning Literature CLS-AN-144 REV2.
55. Rajagopalan LE, and Malter JS, et al. (2000) *J Neurochem.* 74:52-9.
56. Zaric J, and Rüegg C, et al. (2005) *J Biol Chem.* 280:1077-85.
57. Grant DS, et al. (1991) *In Vitro Cell Dev Biol.* 27A:327-36.
58. Reinmuth N, et al. (2003) *Cancer Res.* 63:2079-87.
59. Moser R, et al. (1992) *Blood.* 79:2937-45.
60. Conley BA, et al. (2004) *J. Biol. Chem.* 26:27440-49.
61. Steele JG, et al. (1995) *Biomaterials.* 16:1057-67.
62. Merhi-Soussi F, et al. (2000) *J Leukoc Biol.* 68:881-9.
63. Tanaka Y, et al. (2007) *Clin Cancer Res.* 13:816-23.
64. Panorchan P, et al. (2006) *Biophys J.* 91:3499-507.
65. Tyagi SC. (1998) *Am J Physiol.* 274:C396-405.
66. Saarialho-Kere UK, et al. (1993) *J Clin Invest.* 92:2858-66.
67. Honda K, et al. (2007) *Carcinogenesis.* 28:2074-81.
68. Huang X, et al. (1998) *J Cell Sci.* 111:2189-95.
69. Ito A, et al. (2004) *Tissue Eng.* 10:873-80.
70. Ingthorsson S, et al. (2010) *BMC Res Notes.* 3:184.
71. Miroshnikova YA, et al. (2011) *Phys Biol.* 8:026013.
72. Barcellos-Hoff MH, et al. (1989) *Development.* 105:223-35.
73. Guzman RC, et al. (1987) *Cancer Res.* 47:275-80.
74. Li N, et al. (2005) *EMBO J.* 24:1942-53.
75. Dorosko SM, et al. (2012) *J Pharm Sci.* 101:4666-8.
76. Minakawa T. (1989) *Stroke.* 20:947-51.
77. Tilling T, et al. (1998) *J Neurochem.* 71:1151-7.
78. Gargett CE, et al. (2000) *Hum Reprod.* 15:293-301.
79. Bowman PD, et al. (1982) *In Vitro.* 18:626-32.
80. Belvisi L, et al. (2005) *Mol Cancer Ther.* 4:1670-80.
81. Guo D, et al. (1995) *J Biol Chem.* 270:6729-33.
82. Narmoneva DA, et al. (2005) *Biomaterials.* 26:4837-46.
83. Choi J, et al. (1995) *J Biomed Sci.* 2:379-83.
84. Graeler M, et al. (2002) *J Immunol.* 169:4084-7.
85. Hartley RS, and Yablonka-Reuveni Z. (1990) *In Vitro Cell Dev Biol.* 26:955-61.
86. Kühl U, et al. (1986) *Dev Biol.* 117:628-35.
87. Kase ET, et al. (2007) *Diabetologia.* 50:2171-80.
88. Lamoureux P, et al. (1992) *J Cell Biol.* 118:655-61.
89. Baron-Van Evercooren A, et al. (1982) *J Neurosci Res.* 8:179-93.
90. Armstrong RC, et al. (1997) *J Neurosci.* 17:553 -62.
91. Segal JA, and Skolnick P. (2000) *J Neurochem.* 74:60-9.
92. Wood MW, et al. (2000) *J Neurochem.* 74:2033-40.
93. Malagelada C, et al. (2004) *Stroke.* 35:2396-401.
94. Kushnareva YE, et al. (2005) *J Biol Chem.* 280:28894-902.
95. Lee MK, et al. (1996) *J Neurosci.* 16:7513-25.
96. Tocher DR, and Wilson R. (1990) *J Neurosci Methods.* 33:93-100.
97. Benjamin ER, et al. (2004) *J Biomol Screen.* 9:343-53.
98. Vincent B, et al. (1996) *J Neurosci.* 16:5049-59.
99. Siebler M, et al. (1993) *Brain Res Dev Brain Res.* 73:289-92.
100. Abraham EJ, et al. (2009) Technical Bulletin 471. Corning Literature CLS-DL-CC-004.
101. Laursen LS, et al. (2011) *J Cell Biol.* 192:797-811.
102. Milner R, and French-Constant C. (1994) *Development.* 120:3497-506.
103. Cohen RI, et al. (1996) *J Neurosci.* 16:6433-42.
104. Wozniak M, et al. (2000) *J Bone Miner Res.* 15:1731-45.
105. Bokhari MA, et al. (2005) *Biomaterials.* 26:5198-208.
106. Hulinsky I, et al. (1995) *Pancreas.* 11:309-14.
107. Bonner-Weir S, et al. (2000) *Proc Natl Acad Sci USA.* 97:7999-8004.
108. Buitinga M, et al. (2013) *PLoS One.* 8:e64772.
109. Shitara A, et al. (2009) *Am J Physiol Gastrointest Liver Physiol.* 297:G1198-205.
110. Yeh C, et al. (1991) *In Vitro Cell Dev Biol.* 27A:707-12.
111. Kim MS, et al. (2005) *J Biol Chem.* 280:16163-9.
112. Colombo E, et al. (2013) *PLoS One.* 8:e66734.


References (continued)

113. Pozzetto U, et al. (2005) *Immunol Cell Biol.* 83:25-32.
114. Mavromatis K, et al. (2012) *J Biomol Screen.* 17:1128-35.
115. Mauch S, et al. (2002) *J Immunol.* 168:1244-51.
116. Park JK, et al. (2012) *PLoS One.* 7:e42659.
117. Risner, JR and Holt JR. (2006) *J Neurophysiol.* 96:2364-76.
118. Farrer RG, and Quarles RH. (1996) *J Neurosci Res.* 45:248-57.
119. Moradi F, et al. (2012) *J Neurosci Res.* 90:2335-48.
120. Arnoult C, et al. (1996) *Proc Natl Acad Sci USA.* 93:13004-9.
121. Dym M, et al. (1991) *Endocrinology.* 128:1167-76.
122. Mruk DD, and Cheng CY. (2011) *Methods Mol Biol.* 763:237-52.
123. Funanage VL, et al. (1992) *J Cell Physiol.* 150:251-7.
124. Doi K, et al. (2004) *Interact Cardiovasc Thorac Surg.* 3:359-62.
125. Wensaas AJ, et al. (2007) *J Lipid Res.* 48:961-7.
126. Clyman RI, et al. (1994) *J Cell Sci.* 107:1007-18.
127. Absher M. and Baldor L, et al. (1991) *In Vitro Cell Dev Biol.* 27:25-30.
128. Hedin U, et al. (1988) *J Cell Biol.* 107:307-19.
129. Rieusset J, et al. (2004) *Diabetes.* 53:2232-41.
130. Hudson AE, et al. (2007) *Tissue Eng.* 13:2219-25.
131. Solomon LZ, et al. (1998) *J Lab Clin Med.* 132:279-83.
132. Franck D, et al. (2013) *PLoS One.* 8:e56237.
133. Masters KS, et al. (2004) *J Biomed Mater Res A.* 71:172-80.
134. Rajala K, et al. (2007) *Hum Reprod.* 22:1231-8.
135. Melkounian Z, et al. (2010) *Nat Biotechnol.* 28:606-10.
136. Schulz TC, et al. (2012) *PLoS One.* 7:e37004.
137. Narkilahti S, et al. (2007) *Biomed Eng Online.* 6:11.
138. Jin S, et al. (2012) *PLoS One.* 7:e50880.
139. Bosnakovski D, et al. (2012) Application Note 490. Corning Literature CLS-DL-CC-100.
140. Schwab EH, et al. (2013) *Int J Med Sci.* 10:1846-59.
141. Thompson K, et al. (2012) Application Note 494. Corning Literature CLS-DL-CC-009.
142. Dolley-Sonneville PJ, et al. (2013) *PLoS One.* 8:e70263.
143. Baranov P, et al. (2013) *Tissue Eng Part C Methods.* 19:265-70.
144. Bahrami SB, et al. (2012) *Methods Mol Biol.* 916:81-96.
145. Abe Y, et al. (2013) *PLoS One.* 8:e82085.
146. Marrotte EJ, et al. (2010) *J Clin Invest.* 120:4207-19.
147. Belkind-Gerson J, et al. (2013) *Neurogastroenterol Motil.* 25:61-9.e7.
148. Lancaster MA, et al. (2013) *Nature.* 501:373-9.
149. Yao S, et al. (2006) *Proc Natl Acad Sci USA.* 103:6907-12.
150. Shirasawa S, et al. (2011) *Stem Cells Dev.* 20:1071-8.
151. Laflamme MA, et al. (2007) *Nat Biotechnol.* 25:1015-24.
152. Carpenter L, et al. (2012) *Stem Cells Dev.* 21:977-86.
153. Zhang J, et al. (2012) *Circ Res.* 111:1125-36.
154. Kattman SJ, et al. (2011) *Cell Stem Cell.* 8:228-40.
155. Turovets N, et al. (2012) *Cell Transplant.* 21:217-34.
156. Tada S, et al. (2005) *Development.* 132:4363-74.
157. Rowland TJ, et al. (2010) *Stem Cells Dev.* 19:1231-40.
158. Kopper O, and Benvenisty N. (2012) *Stem Cell Res.* 8:335-45.
159. Lian X, et al. (2013) *Nat Protoc.* 8:162-75.
160. Chen HF, et al. (2007) *Hum Reprod.* 22:567-77.
161. Lyashenko N, et al. (2011) *Nat Cell Biol.* 13:753-61.
162. Fujiwara H, et al. (2007) *J Biol Chem.* 282:29701-11.
163. Tucker BA, et al. (2013) *Stem Cells Transl Med.* 2:16-24.
164. Chiang PM, and Wong PC. (2011) *Development.* 138:2833-43.
165. Battista S, et al. (2005) *Biomaterials.* 26:6194-207.
166. Philp D, et al. (2005) *Stem Cells.* 23:288-96.
167. Joo HJ, et al. (2011) *Blood.* 118:2094-104.
168. Spence JR, et al. (2011) *Nature.* 470:105-9.
169. Sato T, et al. (2011) *Gastroenterology.* 141:1762-72.
170. Carlson AL, et al. (2012) *FASEB J.* 26:3240-51.
171. Chambers SM, et al. (2009) *Nat Biotechnol.* 27:275-80.
172. Hu BY, et al. (2010) *Proc Natl Acad Sci USA.* 107:4335-40.
173. Itsykson P, et al. (2005) *Mol Cell Neurosci.* 30:24-36.
174. Di Giorgio FP, et al. (2008) *Cell Stem Cell.* 3:637-48.
175. Ylä-Outinen L, et al. (2012) *J Tissue Eng Regen Med.* <http://dx.doi.org/10.1002/term.1512>
176. Kim H, et al. (2012) *Stem Cells.* 30:2820-9.
177. Lee TJ, et al. (2013) *Biochem Biophys Res Commun.* 430:793-7.
178. Park SW, et al. (2010) *Blood.* 116:5762-72.
179. Jensen T, et al. (2012) *Tissue Eng Part C Methods.* 18:632-46.
180. Wang D, et al. (2007) *Proc Natl Acad Sci USA.* 104:4449-54.
181. Banerjee ER, et al. (2012) *PLoS One.* 7:e33165.
182. Longmire TA, et al. (2012) *Cell Stem Cell.* 10:398-411.
183. Ishii T, et al. (2008) *Am J Physiol Gastrointest Liver Physiol.* 295:G313-21.
184. Roelandt P, et al. (2010) *PLoS One.* 5:e12101.
185. Sullivan GJ, et al. (2010) *Hepatology.* 51:329-35.
186. Si-Tayeb K, et al. (2010) *Hepatology.* 51:297-305.
187. Yamada T, et al. (2002) *Stem Cells.* 20:146-54.
188. Uemura M, et al. (2010) *J Neurosci Res.* 88:542-51.
189. Ortinau S, et al. (2010) *Biomed Eng Online.* 9:70.
190. Liedmann A, et al. (2012) *Biores Open Access.* 1:16-24.
191. Xia Y, et al. (2013) *Nat Cell Biol.* 15:1507-15.
192. Morizane R, et al. (2013) *PLoS One.* 8:e64843.
193. Guo XM, et al. (2006) *Circulation.* 113:2229-37.
194. Miskon A, et al. (2009) *J Artif Organs.* 12:111-7.
195. Ao A, et al. (2012) *PLoS One.* 7:e41627.
196. Levenberg S, et al. (2003) *Proc Natl Acad Sci USA.* 100:12741-6.
197. Gaspard N, et al. (2008) *Nature.* 455:351-57.
198. Ishizuka T, et al. (2012) *Neurosci Lett.* 525:60-5.
199. Li R, and Mather JP. (2008) *Methods Cell Biol.* 86:227-40.
200. Hutton SR, and Pevny LH. (2008) *CSH Protoc.* Nov 1;2008:pdb.prot5077.
201. Koehler KR, et al. (2013) *Nature.* 500:217-21.
202. Mihaila SM, et al. (2013) *Tissue Eng Part A.* 19:235-46.
203. Wang YH, et al. (2013) *Cytotherapy.* Oct 9. pii: S1465-3249(13)00644-0. doi: 10.1016/j.jcyt.2013.07.015. [Epub ahead of print].
204. Messina E, et al. (2004) *Circ Res.* 95:911-21.
205. Heydarkhan-Hagvall S, et al. (2012) *Biomaterials.* 33:2032-40.
206. Bauer M, et al. (2012) *PLoS One.* 7:e50491.
207. Yui S, et al. (2012) *Nat Med.* 18:618-23.
208. Jung P, et al. (2011) *Nat Med.* 17:1225-27.
209. Vinci M, et al. (2012) *BMC Biol.* 10:29.
210. Yu H, et al. (2006) *Am J Pathol.* 168:1879-88.
211. Lucendo-Villarin B, et al. (2012) *Biores Open Access.* 1:50-3.
212. Soto-Gutiérrez A, et al. (2007) *Nat Protoc.* 2:347-56.
213. Ootani A, et al. (2009) *Nat Med.* 15:701-6.
214. Sato T, et al. (2009) *Nature.* 459:262-5.
215. VanDussen KL, et al. (2012) *Development.* 139:488-97.
216. Snippert HJ, et al. (2012) *Cell.* 143:134-44.
217. Papini S, et al. (2003) *Stem Cells.* 21:481-94.
218. Kim CF, et al. (2005) *Cell.* 121:823-35.
219. Murakami A, et al. (2014) *PLoS One.* 9:e86459.
220. Dong Q, et al. (2013) *Stem Cell Res.* 10:396-404.
221. Guo W, et al. (2012) *Cell.* 148:1015-28.
222. Dey D, et al. (2009) *PLoS One.* 4:e5329.
223. Liu J, et al. (2012) *Stem Cells Dev.* 21:2095-110.
224. Lee JH, et al. (2011) *J R Soc Interface.* 8:998-1010.
225. Qian L. and Saltzman WM. (2004) *Biomaterials.* 25:1331-7.
226. Du J, et al. (2011) *Proc Natl Acad Sci USA.* 108:9466-71.
227. Engler AJ, et al. (2006) *Cell.* 126:677-89.
228. Zhu MX, et al. (2014) *Zhongguo Shi Yan Xue Ye Xue Za Zhi.* 22:148-53.
229. Jiang Y, et al. (2002) *Nature.* 418:41-9.
230. Flanagan LA, et al. (2002) *Neuroreport.* 13:2411-5.
231. Singh P, and Schwarzbauer JE. (2012) *J Cell Sci.* 125:3703-12.

232. Kundu AK, and Putnam AJ. (2006) *Biochem Biophys Res Commun.* 347:347-57.
233. Erickson IE, et al. (2009) *Tissue Eng Part A.* 15:1041-52.
234. Hamada K, et al. (2008) *J Biomed Mater Res A.* 84:128-36.
235. Ozeki M, et al. (2011) *J Biomater Appl.* 25:663-84.
236. Sharma MB, et al. (2012) *Haematologica.* 97:651-60.
237. Chen J, et al. (2013) *Tissue Eng Part A.* 19:716-28.
238. Gutierrez GM, et al. (2008) *Proc Natl Acad Sci USA.* 105:18402-7.
239. Santos Fd, et al. (2011) *Tissue Eng Part C Methods.* 17:1201-10.
240. Allen P, et al. (2011) *J Tissue Eng Regen Med.* 5:e74-86.
241. Gilbert PM, et al. (2010) *Science.* 329:1078-81.
242. Thonhoff JR, et al. (2008) *Brain Res.* 1187:42-51.
243. Addington CP, et al. (2014) *Biomaterials.* 35:3263-72.
244. Lööv C, et al. (2012) *PLoS One.* 7:e29771.
245. Liedmann A, et al. (2012) *J Vis Exp.* Jan 11;(59):e3830.
246. Gelain F, et al. (2006) *PLoS One.* 1:e119.
247. Cunha C, et al. (2011) *Int J Nanomedicine.* 6:943-55.
248. Gao R, et al. (2003) *Diabetes.* 52:2007-15.
249. Yebra M, et al. (2011) *PLoS One.* 6:e22750.
250. Seaberg RM, et al. (2004) *Nat Biotechnol.* 22:1115-24.
251. Georges PC, et al. (2006) *Biophys J.* 90:3012-8.
252. Phillips MJ, and Otteson DC. (2011) *Invest Ophthalmol Vis Sci.* 52:1439-49.
253. Okumura K, et al. (2012) *Stem Cells Int.* 2012:502136.
254. Barker N, et al. (2010) *Cell Stem Cell.* 6:25-36.
255. Uchiyama H, et al. (1992) *Blood.* 80:2306-14.
256. Kessler O, et al. (2004) *Cancer Res.* 64:1008-15.
257. Thompson KE, et al. (2010) Application Note 483. Corning Literature CLS-DL-CC-005.
258. Abraham EJ, and Thompson K. (2009) Application Note 470. Corning Literature CLS-DL-CC-003.
259. Theocharous P, et al. (1997) *J Hematother.* 6:21-9.
260. Schedlich LJ, et al. (1998) *J Biol Chem.* 273:18347-52.
261. Pandrangi SL, et al. (2014) *Cancer Cell Int.* 14:14.
262. Chimal-Ramirez GK, et al. (2013) *Biomed Res Int.* 2013:279505.
263. Ogilvie M, et al. (2000) *J. Biol. Chem.* 275:39754-61.
264. Balci B, and Dinçer P. (2009) *Biotechnol J.* 4:1042-5.
265. Molnar P, et al. (2007) *Biotechnol Prog.* 23:265-8.
266. Yamamoto Y, et al. (2009) *J Biosci Bioeng.* 108:538-43.
267. Kurita-Ochiai T, et al. (2004) *Infect Immun.* 72:5947-54.
268. Servetnyk Z, et al. (2006) *Respir Res.* 7:124.
269. Nishimura S, et al. (2008) *J Biol Chem.* 283:11752-62.
270. Egger M, et al. (1999) *Am J Physiol.* 276:C1184-92.
271. Song H, et al. (2005) *J Virol.* 79:10289-99.
272. Hirahara N, et al. (2000) *Anticancer Res.* 20:3137-42.
273. He J, and Baum LG. (2004) *J Biol Chem.* 279:4705-12.
274. Perry SJ, et al. (2005) *J Biol Chem.* 280:11560-8.
275. Nandivada H, et al. (2012) Application Note 495. Corning Literature CLS-DL-CC-010.
276. Gu M, et al. (1996) *Infect Immun.* 64:524-7.
277. Kim HY, et al. (2007) *J Investig Dermatol Symp Proc.* 12:48-9.
278. Camaire J, and Upton T. (2007) Application Note 1.0.
279. Yano H, et al. (2000) *Proc Natl Acad Sci USA.* 97:9076-81.
280. Thompson EJ, et al. (2006) *Mol Cancer Res.* 4:489-97.
281. Kirsch KH, et al. (1999) *Proc Natl Acad Sci USA.* 96:6211-6.
282. Hu F, et al. (2005) *J Neurosci.* 25:5298-304.
283. Wang S, et al. (1996) *J Biol Chem.* 271:14124-33.
284. Zeng G, et al. (2000) *Cancer Res.* 60:6670-6.
285. Hajj GN, et al. (2007) *J Cell Sci.* 120:1915-26.
286. Valle CW, et al. (2011) *PLoS One.* 6:e29073.
287. Lavoie JN, et al. (2000) *J Cell Biol.* 150:1037-56.
288. Berwin B, et al. (2003) *EMBO J.* 22: 6127-36.
289. Liu J, et al. (2003) *Mol Endocrinol.* 17:346-55.
290. Smith JS, et al. (2001) *J Neurosci.* 21:1096-103.
291. Ho DH, et al. (2004) *Cancer Res.* 64:8968-72.
292. Qian Y, et al. (1997) *J Neurosci.* 17:45-57.
293. Timm DM, et al. (2013) *Sci Rep.* 3:3000.
294. Hearn MG, et al. (2002) *Proc Natl Acad Sci USA.* 99:14554-9.
295. Tom R, et al. (2008) *CSH Protoc.* Mar 1;2008:pdb.prot4978.
296. Pardo AMP, et al. (2010) Technical Report. Corning Literature CLS-AN-057 REV3.
297. Abraham EJ, et al. (2009) Application Note 466. Corning Literature CLS-DL-CC-002.
298. Tian Y, et al. (2013) *PLoS One.* 8:e84457.
299. Corstorphine L, and Sefton MV. (2011) *J Tissue Eng Regen Med.* 5:119-29.
300. Li YW, et al. (2014) *PLoS One.* 9:e87505.
301. Ghisolfi L, et al. (2012) *PLoS One.* 7:e43628.
302. Huyck L, et al. (2012) *Assay Drug Dev Technol.* 10:382-92.
303. Kwak HJ, et al. (2006) *Mol Cancer Res.* 4:209-20.
304. Ishida K, et al. (2004) *Chem Biol.* 11:367-77.
305. Vandenberg P, et al. (1991) *J Cell Biol.* 113:1475-83.
306. Singh S, et al. (2010) *Mol Cancer Ther.* 9:3330-41.
307. Evensen NA, et al. (2013) *PLoS One.* 8:e82811.
308. Liu WD, et al. (2012) *Oncol Lett.* 4:1315-1320.
309. Subramaniam M, et al. (2002) *J Cell Biochem.* 87:9-15.
310. Hosono K, et al. (2007) *Am J Pathol.* 171:274-86.
311. Cuddihy MJ, et al. (2013) *Small.* 9:1008-15.
312. Rivera-Chacon DM, et al. (2013) *J Biomed Nanotechnol.* 9:1092-7.
313. Jang JH, and Kim JH. (2005) *Biotechnol Lett.* 27:1767-70.
314. Zhang B, et al. (2013) *Oncol Lett.* 5:911-916.
315. Thompson K, et al. (2012) Application Note 493. Corning Literature CLS-DL-CC-008.
316. Tschardt M, et al. (2005) *FASEB J.* 19:1836-8.
317. Tiwari A, et al. (2013) *Biopolymers.* 99:334-41.
318. Cheng D, et al. (2008) *J Biomed Mater Res A.* 87:321-31.
319. Vallés-Lluch A, et al. (2013) *Acta Biomater.* 9:9451-60.
320. Zou M, et al. (2013) *Oncol Rep.* 30:615-22.
321. Lee MY, et al. (2008) *Proc Natl Acad Sci USA.* 105:59-63.
322. Doerr ME, and Jones JI. (1996) *J Biol Chem.* 271:2443-7.
323. Komatsu M, et al. (1997) *J Biol Chem.* 272:33245-54.
324. Ferrand N, et al. (2014) *PLoS One.* 9:e87878.
325. Lance A, et al. (2013) *J Tissue Eng Regen Med.* Feb 12. doi: 10.1002/term.1675. [Epub ahead of print].
326. Ampuja M, et al. (2013) *BMC Cancer.* 13:429.
327. Zhang Y, et al. (2011) *J Biol Chem.* 286:16218-28.
328. Adler JJ, et al. (2013) *Proc Natl Acad Sci USA.* 110:17368-73.
329. Miller KA, et al. (2000) *J Biol Chem.* 275:8176-82.
330. Fournier AK, et al. (2008) *J Cell Sci.* 121:226-33.
331. Salas PJ, et al. (1997) *J Cell Biol.* 137:359-75.
332. Guo M, et al. (2002) *J Biol Chem.* 277:14829-37.
333. Wang S, et al. (2013) *J Biol Chem.* 288:26265-74.
334. Yoneda T, et al. (1997) *J Clin Invest.* 99:2509-17.
335. Poincloux R, et al. (2011) *Proc Natl Acad Sci USA.* 108:1943-8.
336. Barneh F, et al. (2013) *Dis Markers.* 35:933-8.
337. Wu GS, et al. (2013) *PLoS One.* 8:e76620.
338. Muñoz-Nájjar UM, et al. (2006) *Oncogene.* 25:2379-92.
339. Hyunju Kim, et al. (2012) *PLoS One.* 7:e49186.
340. Brisson M, et al. (2004) *Mol Pharmacol.* 66:824-33.
341. Li JT, et al. (2013) *Zhonghua Yi Xue Za Zhi.* 93:89-92.
342. Francescone RA 3rd, et al. (2011) *J Vis Exp.* 55:e3040, DOI: 10.3791/3040. URL: <http://www.jove.com/details.php?id=3040>
343. Felding-Habermann B, et al. (2001) *Proc Natl Acad Sci USA.* 98:1853-8.
344. Fabre-Suver C, and Hauschka SD. (1996) *J Biol Chem.* 271:4646-52.
345. Fridman R, et al. (1992) *Int J Cancer.* 51:740-4.
346. Shaw RJ, et al. (1998) *J Biol Chem.* 273:7757-64.
347. Chambers AF, et al. (1993) *Cancer Res.* 53:701-6.

References (continued)

348. Alvarez CV, et al. (1995) *J Biol Chem.* 270:16271-6.
349. Nitsch RM, et al. (1996) *J Biol Chem.* 271:4188-94.
350. Cordeiro ML, et al. (2004) *Neuropsychopharmacology.* 29:39-44.
351. Jacovina AT, et al. (2001) *J Biol Chem.* 276:49350-8.
352. Wu Y, et al. (2004) *Mol Biol Cell.* 15:2093-104.
353. Thomas D, et al. (1995) *J Biol Chem.* 270:28924-31.
354. Hillion JA, et al. (2005) *J Cereb Blood Flow Metab.* 25:154-62.
355. Riederer BM, et al. (1997) *Proc Natl Acad Sci USA.* 94:741-5.
356. Holgado-Madruga M, et al. (1997) *Proc Natl Acad Sci USA.* 94:12419-24.
357. Hiromura M, et al. (2006) *J Biol Chem.* 281:27753-64.
358. Jos A, et al. (2005) *Aquat Toxicol.* 71:183-92.
359. McKay DB, et al. (1992) *J Clin Invest.* 89:134-42.
360. Crasto,GJ, et al. (2010). Customer Application Note. Corning Literature CLS-AN-165.
361. Kartner N, et al. (2010) *J Biol Chem.* 285:37476-90.
362. Ivankovic-Dikic I, et al. (2000) *Nat Cell Biol.* 2:574-81.
363. Lu HF, et al. (2009) *Neurochem Res.* 34:1575-83.
364. Leventhal PS and Feldman EL. (1996) *J Biol Chem.* 271:5957-60.
365. Meyer A, et al. (2004) *Neoplasia.* 6:332-42.
366. Nakanon J, et al. (1999) *J Investig Dermatol Symp Proc.* 4:173-6.
367. Conner SR, et al. (2003) *J Biol Chem.* 278:34548-54.
368. Van Goethem E, et al. (2010) *J Immunol.* 184:1049-61.
369. Kurata S, et al. (2004) *J Biol Chem.* 279:50069-77.
370. Waltz DA, et al. (1997) *J Clin Invest.* 100:58-67.
371. Liu Z, et al. (2011) *J Biol Chem.* 286:39654-62.
372. Li L, and Lu Y. (2011) *J Cancer.* 2:458-66.
373. Butcher JT, et al. (2006) *Biomaterials.* 27:5252-8.
374. Hwang DS, et al. (2007) *J Biotechnol.* 127:727-35.
375. Hemmings KE, et al. (2012) *Biol Reprod.* 86:165.
376. Guo C, and Piacentini L. (2003) *J Biol Chem.* 278:46699-708.
377. Zhong Wu X, et al. (2004) *Int J Cancer.* 110:504-10.
378. Williams JL. (1999) *Am J Trop Med Hyg.* 60:7-13.
379. Schuster FL. (2002) *Clin Microbiol Rev.* 15:355-64.
380. Chen T, et al. (2007) *Stem Cells.* 25:392-401.
381. Li Z, et al. (2009) *PLoS One.* 4:e8443.
382. Vörsmann H, et al. (2013) *Cell Death Dis.* 4:e719.
383. Greggio C, et al. (2013) *Development.* 140:4452-62.
384. Touboul C, et al. (2013) *J Transl Med.* 11:28.
385. Notter MF. (1988) *Exp Cell Res.* 177:237-46.
386. Kim Y, and Rajagopalan P. (2010) *PLoS One.* 5:e15456.
387. Yamajuku D, et al. (2012) *Sci Rep.* 2:439.
388. Semino CE, et al. (2003) *Differentiation.* 71:262-70.
389. Wang MF, et al. (2011) *Zhongguo Shi Yan Xue Ye Xue Za Zhi.* 19:944-8.
390. Rah B, et al. (2012) *PLoS One.* 7:e44039.
391. Ramer R, and Hinz B. (2008) *J Natl Cancer Inst.* 100:59-69.
392. Rodin S, et al. (2014) *Nat Commun.* 5:3195.
393. Lu HF, et al. (2014) *Biomaterials.* 35:2816-26.
394. Burridge PW, et al. (2014) *Nat Methods.* 11:855-60.
395. Maruotti J, et al. (2013) *Stem Cells Transl Med.* 2:341-54.
396. Zahabi A, et al. (2012) *Stem Cells and Dev.* 21:2262-72.
397. Pardo AP, and Rothenberg ME. (2013, 2014) Application Note. Corning Literature CLS-AN-227 REV1.

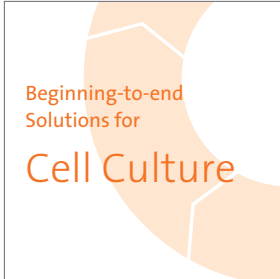

www.corning.com/lifesciences/solutions

At Corning, cells are in our culture. In our continuous efforts to improve efficiencies and develop new tools and technologies for life science researchers, we have scientists working in Corning R&D labs across the globe, doing what you do every day. From seeding starter cultures to expanding cells for assays, our technical experts understand your challenges and your increased need for more reliable cells and cellular material.

It is this expertise, plus a 160-year history of Corning innovation and manufacturing excellence, that puts us in a unique position to offer a beginning-to-end portfolio of high-quality, reliable cell culture consumables.

Corning Incorporated
Life Sciences

836 North St.
Building 300, Suite 3401
Tewksbury, MA 01876
t 800.492.1110
t 978.442.2200
f 978.442.2476
www.corning.com/lifesciences

**Worldwide
Support Offices**

ASIA/PACIFIC
Australia/New Zealand
t 61 427286832
China
t 86 21 3338 4338
f 86 21 3338 4300
India
t 91 124 4604000
f 91 124 4604099

Japan
t 81 3-3586 1996
f 81 3-3586 1291

Korea
t 82 2-796-9500
f 82 2-796-9300

Singapore
t 65 6572-9740
f 65 6861-2913

Taiwan
t 886 2-2716-0338
f 886 2-2516-7500

EUROPE

France
t 0800 916 882
f 0800 918 636

Germany
t 0800 101 1153
f 0800 101 2427

The Netherlands
t 31 20 655 79 28
f 31 20 659 76 73

United Kingdom
t 0800 376 8660
f 0800 279 1117

**All Other European
Countries**
t 31 (0) 20 659 60 51
f 31 (0) 20 659 76 73

LATIN AMERICA
grupoLA@corning.com
Brasil
t (55-11) 3089-7400
f (55-11) 3167-0700
Mexico
t (52-81) 8158-8400
f (52-81) 8313-8589



For a listing of trademarks, visit us at www.corning.com/lifesciences/trademarks.
All other trademarks are the property of their respective owners.