

PYREX[®] Laboratory Glassware

Fuel your expertise with the lab glassware trusted by scientists for more than 100 years.

PYREX[®]

A Corning Brand





Table of Contents

Exceptional strength and stability. Brilliant transparency. Superior resistance to chemicals, contaminants, and drastic changes in temperature. These are among the key properties that enable PYREX® glass to endure even the toughest laboratory environments, making PYREX scientific glassware the indispensable partner for chemists and research scientists for more than 100 years.

PYREX® Glassware	3
PYREX VISTA™ Glassware	4
PYREXPLUS® Glassware	5
PYREX Low Actinic Glassware	6

PYREX[®] Glassware

Of the hundreds of commercial glasses produced, Code No. 7740 borosilicate glass comes closest to being the ideal glass for most laboratory applications. With proper care, it will withstand nearly all temperatures used in normal laboratory use. It is highly resistant to chemical attack. Its low coefficient of expansion allows it to be manufactured to suit a variety of wall thicknesses, tailored to different purposes. This mechanical strength is not solely dependent on wall thickness but is also influenced by the glass formula and type, ensuring durability while maintaining reliable heat resistance.

Ordering Information

Cat. No. (Series)	Volumes	Description	Pack/Case
1000	10 mL - 4L	Beakers, low form, double scale, graduated	6 to 48
1395	25 mL - 10L	Media storage bottles, round, screw cap	1 to 10
3023	25 mL - 2L	Cylinders, Class A, double metric scale, to deliver	1 to 24
4980	10 mL - 6L	Erlenmeyer flasks, narrow mouth, heavy duty rim, graduated	1 to 48
5640	1 mL - 6L	Volumetric flasks, Class A, PYREX standard taper stopper	1 to 12
9820	0.5 mL - 70 mL	Culture tubes, rimless	192 to 720
9825	9 mL - 170 mL	Culture tubes, screw cap	48 to 288



PYREX[®] VISTA[™] Glassware

PYREX VISTA glassware is an economical option for the customer who is willing to forgo the premium benefits of PYREX products. Manufactured to Corning/PYREX standards and priced competitively against comparable products, PYREX VISTA glassware offers a full range of products from beakers to pipets and is easily recognized by its blue graduations and novel marking spot.

Ordering Information

Cat. No. (Series)	Volumes	Description	Pack/Case
70000	10 mL - 1L	Beakers, standard low griffin	12 to 24
70024	10 mL - 500 mL	Cylinders, Class A, single metric scale, graduated, funnel top, to deliver	8 to 24
70980	125 mL - 1L	Erlenmeyer flasks, narrow mouth	12 to 24
70640	10 mL - 2L	Volumetric flasks, Class A, glass standard taper stopper	2 to 12
70800	9 mL - 70 mL	Test tubes, beaded rim	200
70825	9 mL - 55 mL	Culture tubes, screw cap, phenolic cap	200




PYREXPLUS® Glassware

PYREXPLUS glassware is coated with a tough, transparent plastic vinyl. The coating, which is applied to the outside of the vessel, helps prevent exterior surface abrasion. It also helps minimize the loss of contents and helps contain glass fragments if the glass vessel is broken.

Ordering Information

Cat. No. (Series)	Volumes	Description	Pack/Case
61626	100 mL - 2L	Media storage bottles, PVC protective coating, round, screw cap	4
63024	100 mL - 2L	Cylinders, single metric scale, protective coating, graduated, to deliver	2 to 4
65340	250 mL - 4L	Filtering flasks, heavy wall, protective coating, tubulation, graduated	2 to 6
65640	100 mL - 2L	Volumetric flasks, Class A, protective coating, PYREX standard taper stopper	2 to 12
66402	125 mL - 2L	Separatory funnels, Squibb pear-shaped, protective coating, PTFE standard stopcock, glass standard taper stopper	1 to 4

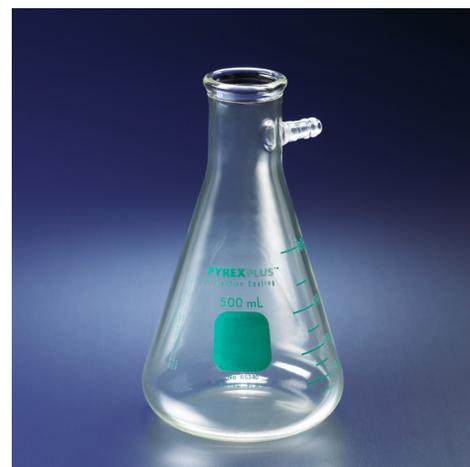


PYREX[®] Low Actinic Glassware

PYREX Low Actinic Labware is made from Code. No. 7740 glass with copper replacing the sodium in the glass matrix at the surface. The resulting product is as durable as the base glass. PYREX Low Actinic Labware was originally developed for work in the vitamin field, but it has found other uses in applications with chemicals sensitive to light in the 3,000 to 5,000 Angstrom range. Low actinic glassware is designed to reduce light transmission, especially UV and visible wavelengths, protecting sensitive samples such as chemicals, reagents, and biological materials from degradation caused by light exposure. By minimizing light-induced reactions, the glassware ensure the stability and integrity of stored or processed substances, which is critical for precision in experiments and long-term storage. PYREX Low Actinic glassware is ideal for use in chemistry, biology, pharmaceuticals, and photochemistry experiments.

Ordering Information

Cat. No. (Series)	Volumes	Description	Pack/Case
51395	25 mL - 10L	Media storage bottles, low actinic, round, graduated, with screw cap	1 to 10
3042	100 mL - 2L	Cylinders, Lifetime Red™ band marking, single metric scale, graduated, to contain	1 to 24
55640	10 mL - 2L	Volumetric flasks, low actinic, Class A, PYREX standard taper stopper	1 to 12
5660	25 mL - 2L	Volumetric flasks, Lifetime Red band marking, Class A, PYREX standard taper stopper, to contain	4 to 12

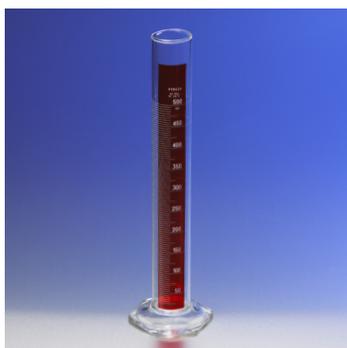

How to Purchase: For specific availability in your region and purchasing options, terms and conditions of sale, customer/product support, and certificates, visit www.corning.com/how-to-buy.

Warranty/Disclaimer: Unless otherwise specified, all products are for research use or general laboratory use only.* Not intended for use in diagnostic or therapeutic procedures. Not for use in humans. These products are not intended to mitigate the presence of microorganisms on surfaces or in the environment, where such organisms can be deleterious to humans or the environment. Corning Life Sciences makes no claims regarding the performance of these products for clinical or diagnostic applications. *For a listing of US medical devices, regulatory classifications or specific information on claims, visit www.corning.com/resources.

Corning's products are not specifically designed and tested for diagnostic testing. Many Corning products, though not specific for diagnostic testing, can be used in the workflow and preparation of the test at the customers discretion. Customers may use these products to support their claims. We cannot make any claims or statements that our products are approved for diagnostic testing either directly or indirectly. The customer is responsible for any testing, validation, and/or regulatory submissions that may be required to support the safety and efficacy of their intended application.

CORNING

Corning Incorporated
Life Sciences

www.corning.com/lifesciences

NORTH AMERICA

t 800.492.1110

t 978.442.2200

ASIA/PACIFIC

Australia/New Zealand

t 61 427286832

Chinese Mainland

t 86 21 3338 4338

India

t 91 124 4604000

Japan

t 81 3-3586 1996

Korea

t 82 2-796-9500

Singapore

t 65 6572-9740

Taiwan

t 886 2-2716-0338

EUROPE

CSEurope@corning.com

France

t 0800 916 882

Germany

t 0800 101 1153

The Netherlands

t 020 655 79 28

United Kingdom

t 0800 376 8660

All Other European Countries

t +31 (0) 206 59 60 51

LATIN AMERICA

grupolA@corning.com

Brazil

t 55 (11) 3089-7400

Mexico

t (52-81) 8158-8400

The information contained within is accurate as of the date of publication and subject to change without notice.

For a listing of trademarks, visit www.corning.com/trademarks. All other trademarks are the property of their respective owners.

© 2025, 2026 Corning Incorporated. All rights reserved. 1/26 CLS-GL-064 REV1